


pCMV-HA-C

Catalog No(s).
635690

pCMV-HA-C Vector Information

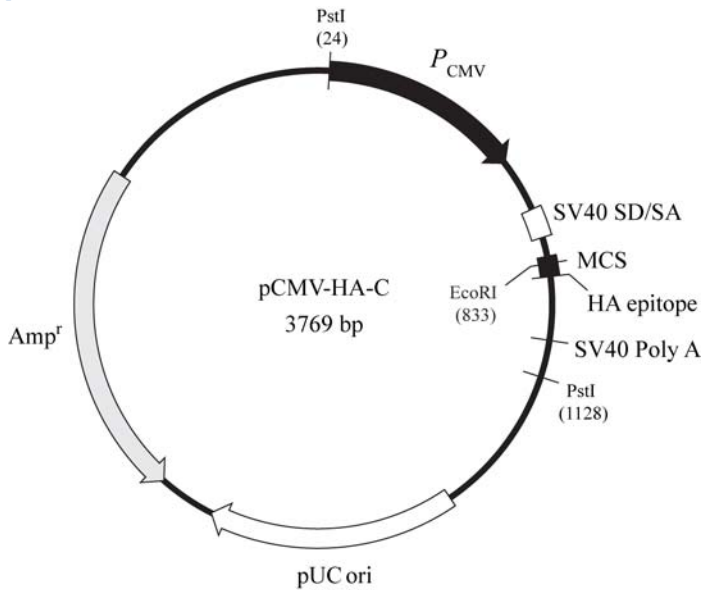

	ApaI		EcoRI		SalI		BglII		XhoI
819	CGG GCC CAG GCC CGA ATT CGG TCG ACC GAG ATC TCT CGA		GCC CGG GTC CGG GCT TAA GCC AGC TGG CTC TAG AGA GCT						
	Acc65I								
	XhoI								
	KpnI								
858	GGT ACC GGC TAC CCA TAC GAT GTT CCA GAT TAC GCT TAA		TAC CCA TAC GAT GTT CCA GAT TAC GCT TAA						
	CCA TGG CCG ATG GGT ATG CTA CAA GGT CTA ATG CGA ATT								

pCMV-HA-C vector map and multiple cloning site.

To create a fusion of your protein of interest and the HA tag, insert your gene of interest in-frame with the HA coding sequence.

Location of Features

- P_{CMV} (human cytomegalovirus promoter): 27–549
 - SV40 SD/SA (SV40 donor/acceptor splice sites): 672–768
 - MCS (multiple cloning site): 816–863
 - HA epitope: 867–896
 - SV40 polyA (SV40 polyA signal): 1031–1050
 - pUC ori (pUC origin of replication): 1523–2166
 - Amp^r (ampicillin resistance gene): 2314–3174 (complementary)
-

Notice to Purchaser

Clontech products are to be used for research purposes only. They may not be used for any other purpose, including, but not limited to, use in drugs, *in vitro* diagnostic purposes, therapeutics, or in humans. Clontech products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without prior written approval of Clontech Laboratories, Inc.

Your use of this product is also subject to compliance with the licensing requirements described on the product's web page at <http://www.clontech.com>. It is your responsibility to review, understand and adhere to any restrictions imposed by these statements.

Clontech and the Clontech logo are trademarks of Clontech Laboratories, Inc. All other marks are the property of their respective owners. Certain trademarks may not be registered in all jurisdictions. Clontech is a Takara Bio Company. ©2011 Clontech Laboratories, Inc.

This document has been reviewed and approved by the Clontech Quality Assurance Department.